



BlueTip.

Application Instructions
Step by Step.



[01] Screw the BlueTip activation tip into the thread of the sonic handpiece by hand in a clockwise direction until you feel a slight increase in resistance. Take care not to

- bend or
- kink the activation tip.



[02] Fill the pulp chamber with one of the following irrigation solutions:

NaOCl solution	max. 5.25%
EDTA	max. 17%
Citric acid	max. 10%
H ₂ O ₂	max. 3%
CHX	max. 2%



[03] Introduce the activation tip into the root canal and switch on the sonic handpiece and leave it on for up to 30 seconds as required, moving the activation tip up and down without applying pressure.

Important:


The frequency of the sonic handpiece must be set to the highest power level, making sure not to exceed the upper limit of 6000 Hz

(Air supply at a maximum of 0.3 MPa (3 bar)).


Only activate the irrigation tip when it is inside the root canal.

The irrigation tip is a sterile, disposable product. It is intended for use on a single patient only. The irrigation tip cannot be disinfected or autoclaved.


Reusing the irrigation tip increases the risk of spontaneous fracture of the tip and may lead to cross-contamination.

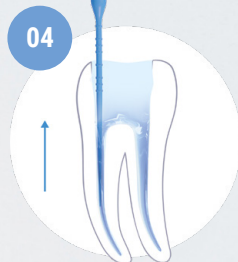

[04] Repeat the irrigation process, activating the irrigation solution with the activation tip after each preparation step.

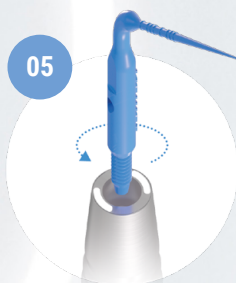

[05] Upon completion of the application: Unscrew the activation tip from the sonic handpiece by hand, turning it in an anti-clockwise direction.

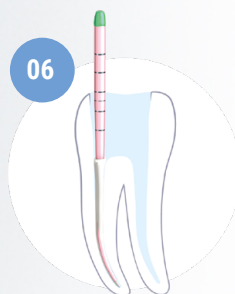


[06] Dry the canal and obturate it using your preferred method, for example the single-cone technique in combination with Komet **BioSeal**.

Only use commercially available irrigation solutions designed for this purpose when irrigating the root canal:

NaOCl solution		max. 5.25%
EDTA		max. 17%
Citric acid		max. 10%
H ₂ O ₂		max. 3%
CHX		max. 2%

Good to know:

compatible sonic handpieces include
Komet SF1LM / SF1LS, W&H ZA-55, NSK
Ti-Max 970, KaVo SONICflex 2003

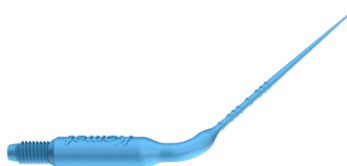

SFPA66.000.020

STERILE



Komet Dental

Gebr. Brasseler GmbH & Co. KG
Trophagener Weg 25 · 32657 Lemgo
Postfach 160 · 32631 Lemgo · Germany
Phone +49 (0) 5261 701-0
hello@kometdental.com

Komet Germany

Trophagener Weg 25 · 32657 Lemgo
Phone +49 (0) 800 7701-700
Fax +49 (0) 800 7701-800
info@kometdental.de
www.kometdental.de

Komet Austria

Komet Austria Handelsagentur GmbH
Hellbrunner Straße 15 · 5020 Salzburg
Phone +43 (0) 0800 70170070
Fax +43 (0) 662 829-435
info@kometdental.at
www.kometdental.at

Komet France

SA à conseil d'administration, s.a.i.
18 rue · Fourcroy · 75017 Paris
Phone +33 (0) 1 43 48 89 90
info@komet.fr
www.komet.fr

Komet Italy

Komet Italy S.r.l.
Via Gianbattista Morgagni 36,
37 135 Verona
Phone +39 (0) 45 11171911
info@komet.it
www.komet.it

Komet UK

Komet Dental Instruments UK Limited
Unit 8 Arden House · Marsh Lane
Hampton In Arden · B92 0AJ
Phone +44 (0) 2046 200580
info@kometuk.com
www.kometuk.com

Komet International Distribution

Trophagener Weg 25 · 32657 Lemgo
Phone +49 (0) 5261 701-0
export@kometdental.de
www.kometdental.com

