




Komet **EnFill & EnPack.**

Outstanding on their own. Even better as a team.



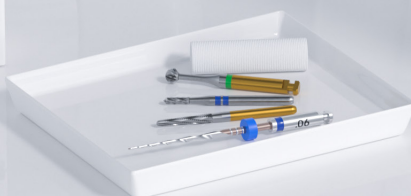

Komet **EnFill & EnPack.**

The combination for economical, simple and safe root canal obturations.

Highly efficient, safe and homogeneous root canal filling - the perfect end of every endodontic treatment. Battery-operated, precise, ready for use in a flash

With our obturation devices for warm, vertical obturation, you have the ideal combination in your hands to complete your endodontic treatments in an economically sound, comfortable and safe manner.

With EnPack and EnFill, you are always ready for action. You can get started straight away thanks to the extremely short heating times and work efficiently without interruptions.

Powerful rechargeable batteries without annoying cables, ease of operation thanks to the plug-and-play principle and the unique, rotatable displays ensure that you can concentrate fully on your work without worrying about the technology. A quick glance, a quick push of a button, a quick heat-up and you're ready to go.

EnPack and EnFill have everything you need to maintain your usual workflow. Both devices can be used universally, they meet almost any requirement thanks to the 360° rotatable injection needle and their heating tips in nine sizes. The devices allow you to work without fatigue thanks to their ergonomic design and low weight.

Increase your treatment efficiency with outstanding, homogeneous results in record time, from the apex to the canal entrance.

Komet **EnPack.**

Efficiency in the root canal.

With our obturation device for warm, vertical obturation of the apical third, you can complete your endodontic treatments economically, comfortably and safely.

You are ready for action at any time, can get started straight away and work efficiently without interruptions. The heating tip heats up in just 0.2 seconds, so you don't lose any valuable time during work and stay in the workflow. Simple plug-and-play operation, a display that is easy to read in any position and a powerful battery ensure that you only have to worry about your patient and not the technology.

EnPack makes clean, precise work particularly easy. It is suitable for universal use, i.e. by those who prefer a purely manual preparation as well as in conjunction with any endodontic motor. The nine heating tips with different diameters always offer an optimum view of the working area.

Thanks to the low weight and ergonomic design of the device, you can work without fatigue. Precise separation of the master point is easy with the heating tip even when working with cold obturations - without clogging up and the risk of pulling the gutta-percha out again.

Rapid heating in the root canal ensures maximum patient safety and a perfect, homogeneous result.

Use:

1. Briefly press the power button to switch on the EnPack.
2. Press and hold the power button to change the display.
3. Press the heating button to heat the heating tip. Release it to end the heating process.
4. After switching on, press the power button and the heating button simultaneously to adjust the volume.
5. Briefly press the power button to select the volume.
6. Briefly press the Heating button to end the volume setting.





Your highlights. Our promise.

Everything at a glance.

Ready at any time.

The heating tip heats up in just 0.2 seconds, so you can get started straight away. The tip heats up directly inside the canal, so that you increase the safety of your patients and ensure optimal processing of the gutta-percha.

Save time with every treatment.

No lengthy preparation, just get started: With just two push buttons and a few settings, you can simply do your job without worrying about the technology. The simple operation allows you to react quickly and save valuable time with every treatment.

Full focus on what really matters.

The tip comes in different diameters and offers an optimum view of the working area at all times. It reaches even narrow, deep areas and separates the gutta-percha with the utmost accuracy. This allows you to work in comfort and with maximum precision, without restrictions even in confined situations.

Equipped for every situation.

With nine heating tip sizes and a precisely adjustable angle, almost any root canal can be filled efficiently. In just a few simple steps, your EnPack is ready for the next patient - precisely tailored to the requirements of the root canal.

Work comfortably.

Work without fatigue: Thanks to its slim design, low weight and ergonomic handle, endodontic treatments with the EnPack can be done in perfect comfort. The design and display are suitable for both left-handed and right-handed users.

Highly secure.

Fast heating time and rapid cooling offer the greatest possible safety for the patient and perfect results. The patient is additionally protected because the tip is heated up inside the root canal.

Versatile compatibility.

As a stand-alone device, the EnPack is the perfect addition to any endo motor. Suitable for universal use, it pairs up ideally with the EnFill and perfectly complements our EnDrive.

Universal use.

You can also benefit from the precision of the EnPack during cold obturation where it separates the master point cleanly and precisely. No clogging up, no risk of moving the gutta-percha, no risk of injury, just comfortable work with a sterile, clean result.

Flexible use.

No annoying cables, just a battery that lasts for 1,500 treatments or 4 hours. When not in use, simply place your EnPack on the charging station and it charges up automatically.

Komet **EnFill.**

Simple and safe filling.

With our obturation device for warm, vertical obturations, you can complete your endodontic treatments in an economically sound, comfortable and safe manner. Thanks to the powerful rechargeable battery and the short heating time of only around 15 seconds, you are ready to go in no time at all and can carry out your filling without interruptions or delays.

The six available injection needles can be optimally adapted to the requirements of the root canal and your preferences, as they can be rotated by 360°. They are so thin that you have a perfect view of the working area at all times and can fill the canal safely and precisely. Intuitive operation, fast responsiveness and fatigue-free design save valuable time during every endodontic treatment. The EnFill can be used by those who prefer manual preparations as well as in combination with any endodontic motor.

The lightweight, ergonomic pistol grip gives you full control over the filling volume and flow speed and follows your own workflow. Whether on its own or in combination with the EnPack as an ideal partner – with EnFill you will also reap the rewards for all the challenging work you have already done with your root canal preparation.

Use:

1. Press and hold the release mechanism and pull back the tappet.
2. Fill the loading slot with a gutta-percha pellet. **IMPORTANT:** Only ever use ONE pellet!
3. Push the tappet back into the appliance until the black mark disappears.
4. Press and hold the power button to switch on the EnFill. .
5. When the appliance is switched on, briefly press the power button to set the desired temperature (recommendation: 180° C or 200° C).
6. Once the desired temperature has been reached, gently press the tappet forward to inject the filling compound.



Why it pays off for you.

Everything at a glance.

Stay efficient in the workflow.

Stay in the workflow instead of waiting for your device. The EnFill heats the gutta-percha to the desired temperature in just 15 seconds. This allows you to continue working almost seamlessly, without interruptions.

Meets your requirements exactly.

Every root canal is unique, just like your habits. That's why each of the 6 available injection needles can be freely rotated by 360°, thus perfectly adapting to the changing requirements of your endodontic treatment.

Full focus on what really matters.

The EnFill injection needles come in different diameters so that you can safely reach almost any area that needs to be filled. They are so thin that you always have the best view of the working area, no matter how tight it gets.

Save time with every treatment.

With just one push button and a few settings, you can simply get on with your work without having to worry about the technology. The easy operation allows you to react quickly and save valuable time with every treatment.



Work comfortably and precisely.

Precise work without interruption or fatigue: The pistol grip is designed in such a way that you can easily and precisely control the flow rate and quantity of gutta-percha. Endodontic treatments can be carried out in complete comfort thanks to the low weight and ergonomic handle of the device.

Particularly safe.

Fast heating-up period, short application time and rapid cooling to ensure the greatest possible safety for the patient and perfect results without loosening of the filling. The patient is additionally protected because the device heats up in the root canal.

Versatile compatibility.

As a stand-alone device, the EnPack is the perfect addition to any endo motor. Suitable for universal use, it pairs up ideally with the EnFill and perfectly complements our EnDrive

Flexible use.

No annoying cables, just a battery that usually lasts for an entire day of use at the practice. When not in use, simply place your EnFill on the charging station and it will charge up automatically.

Komet Dental

Gebr. Brasseler GmbH & Co. KG
Trophagener Weg 25 · 32657 Lemgo
Postfach 160, 32631 Lemgo · Germany
Telefon +49 (0) 5261 701-0
hello@kometdental.com
Contact for dealers:
export@kometdental.de

Verkauf Deutschland:

Telefon +49 (0) 5261 701-700
Telefax +49 (0) 5261 701-289
info@kometdental.de
www.kometdental.de

Export:

Telefon +49 (0) 5261 701-0
Telefax +49 (0) 5261 701-329
export@kometdental.de
www.kometdental.com

Komet Austria Handelsagentur GmbH

Hellbrunner Straße 15 · 5020 Salzburg
Telefon +43 (0) 662 829-434
Telefax +43 (0) 662 829-435
info@kometdental.at
www.kometdental.at

Komet France

18 rue · Fourcroy · 75017 Paris
Phone +33 (0) 1 43 48 89 90
info@komet.fr
www.komet.fr

Komet Italy

Via Gianbattista Morgagni 36,
37135 Verona
Phone +39 (0) 45 11171911
info@komet.it
www.komet.it

Komet UK

Invoice and general correspondence:
27 Old Gloucester St,
Holborn, London WC1N 3AX
info@kometuk.com
www.kometuk.com

Brasseler®, Komet®, Art2®, CeraBur®, CeraCut®, CeraDrill®,
CeraPost®, DC1®, DIAO®, FastFile®, F360®, F6 SkyTaper®,
H4MC®, KometBioSeal®, OccluShaper®, OptiPost®, PolyBur®,
PrepMarker®, Procodile®, Procodile Q® et SHAX® are registered
trademarks of Gebr. Brasseler GmbH & Co. KG.

ReFlex® and EndoPilot® are registered trademarks of
Schlumbohm GmbH & Co. KG

Some of the products and designations mentioned in
the text are trademarked, patented or copyrighted.
The absence of a special reference or the sign ® should not be
interpreted as the absence of legal protection.

This publication is copyrighted. All rights, also with regard to
translation, reprint and reproduction (also in the form of extracts)
are reserved. No part of this publication may be reproduced or
reprocessed using electronic systems in any form or by any means
(photocopying, microfilm or other methods) without the written
permission of the editor.

Colours and products are subject to alterations.
Printing errors excepted.

As at February 2025