




Dentastic Remover **Kit 4687B** - Precision. Speed. Reliability.

Effortless Removal, Exceptional Performance.

Introducing the **Dentastic Remover Kit**, your ultimate solution for precise and efficient removal of a wide range of dental restorations. Designed for professionals who demand performance and reliability, this kit includes a selection of top-tier instruments that ensure smooth, safe, and time-saving procedures. From amalgam fillings to metal and ceramic crowns, the Dentastic Remover Kit equips you with the tools needed to streamline your workflow while ensuring patient comfort.

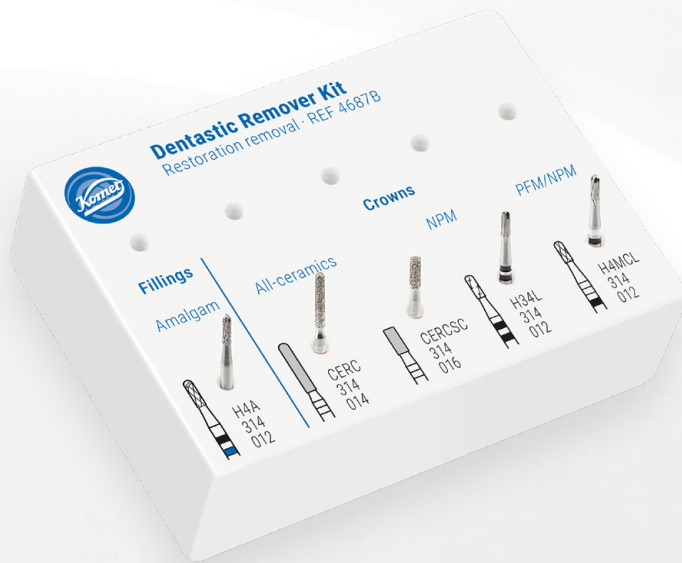



What's Inside?

Content of Kit 4687B:



H4A – Amalgam remover:

Sharper equals faster equals safer.

Increase your treatment efficiency and improve the safety of your patients, your staff and your own: With its new, sharper blade, the H4A amalgam remover reduces the emission of mercury vapors by 40%.*

At the same time, the powerful cutting performance enables controlled, efficient work and protects the natural tooth structure.

H4A – for more protection and efficiency at the dental practice - every day.

*Source: Komet Dental TestLab, mechanical cutting test 2025.
The test is based on the comparison with the predecessor H32.



H34L – The Multi-Talented Metal Cutter

The H34L is our versatile cutter designed for all metal crowns, effortlessly cutting through metal alloys without difficulty.

Featuring precision-engineered cross-cut toothings, it ensures efficient and precise separation while allowing ample space for rapid chip evacuation. Equipped with a pronounced tip-transversing blade, the H34L facilitates smooth penetration of crowns, ensuring seamless operation and long service life.

For optimal performance, the instrument is applied at a 45-degree angle and operated at an ideal speed of 160,000 rpm.

To achieve the best results, work with minimal contact pressure and ensure continuous cooling during separation.



H4MCL – The Crown Removal Specialist

The H4MCL crown remover series stands out with its pyramid-shaped blade configuration, delivering exceptional durability.

Engineered for efficiency, it removes crowns made of metal alloys, titanium, and thin ceramic veneers with low fusion points.

Its unique blade design ensures optimal performance, while the tip-transversing blade allows for easy penetration, helping to streamline procedures and maximise efficiency in your practice.

Available in various working part lengths, the H4MCL is a reliable choice in an ever-changing world.




**Rocky – The Champion in crown-cutting**

Rocky revolutionises all-ceramic crown removal with 74%* greater initial sharpness and 87%* faster cutting.

Even tough zirconia is removed effortlessly, reducing procedure time.

Its innovative next-gen diamond grit ensures exceptional sharpness, durability, and precision for smoother, faster work with less effort. The distinctive rose gold shank allows quick identification in busy practices, blending performance with sophistication.

*Source: Komet Dental TestLab, mechanical cutting test 2025.
Test based on comparison to competitive average

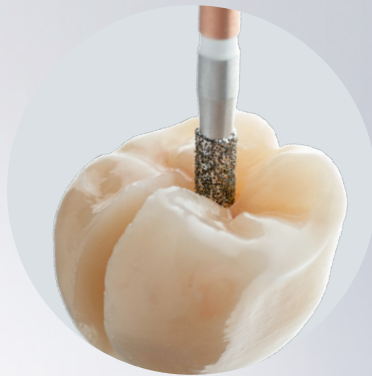

Why Choose the Dentastic Remover Kit?

- **Unmatched Efficiency**
Remove restorations with speed and precision.
- **Enhanced Patient Comfort**
Minimises heat generation and procedure time.
- **Versatile Application**
Handles amalgam, metal, and ceramic restorations.
- **Long-Lasting Performance**
Durable, high-quality instruments for extended use.

Elevate your dental practice with the **Dentastic Remover Kit 4687B** – where innovation meets efficiency.

Komet Dental – Precision Tools. Expert Guidance. Tailored for You.

Komet Dental

Gebr. Brasseler GmbH & Co. KG
Trophagener Weg 25 · 32657 Lemgo
Postfach 160 · 32631 Lemgo · Germany
Phone +49 (0) 5261 701-0
hello@kometdental.com

Komet Germany

Trophagener Weg 25 · 32657 Lemgo
Phone +49 (0) 800 7701-700
Fax +49 (0) 800 7701-800
info@kometdental.de
www.kometdental.de

Komet Austria

Komet Austria Handelsagentur GmbH
Hellbrunner Straße 15 · 5020 Salzburg
Phone +43 (0) 0800 70170070
Fax +43 (0) 662 829-435
info@kometdental.at
www.kometdental.at

Komet France

SA à conseil d'administration, s.a.i.
18 rue · Fourcroy · 75017 Paris
Phone +33 (0) 1 43 48 89 90
info@komet.fr
www.komet.fr

Komet Italy

Komet Italy S.r.l.
Via Gianbattista Morgagni 36,
37 135 Verona
Phone +39 (0) 45 11171911
info@komet.it
www.komet.it

Komet UK

Komet Dental Instruments UK Limited
Unit 8 Arden House · Marsh Lane
Hampton In Arden · B92 0AJ
Phone +44 (0) 2046 200580
info@kometuk.com
www.kometuk.com

Komet International Distribution

Trophagener Weg 25 · 32657 Lemgo
Phone +49 (0) 5261 701-0
export@kometdental.de
www.kometdental.com